

Super Mouse 1800™ Stack-A-Rack™


Maximum Flexibility

Create your own stackable housing system with our unique Stack-A-Rack™ Interchangeable See-Through™ Housing Systems. Designed specifically for housing mice, the Super Mouse 1800™ Stack-A-Rack™ provides the ability to house both the Super Mouse 1800 cage and the smaller Super Mouse 750™ cages, for maximum flexibility. Simply rotate cage runner supports to switch from larger to smaller cages and house both simultaneously on the same rack. Up to eight modular tiers accept as many as 64 Super Mouse 1800™ or 128 Super Mouse 750™ cages—or anywhere in between!

Stack-A-Rack™ systems provide a variety of ordering and configuration options. To meet various housing density and installation requirements, single-sided and double-sided mobile systems are available. You can order Stack-A-Rack™ unassembled, preassembled or with welded, fixed shelving. Super Mouse 1800™ and Super Mouse 750™ cages accept LifeSpan™ rodent enrichment, providing a stimulating multilevel cage environment. Both cages also accept Micro-Isolator™ Filter Tops, allowing cage inventory to be shared with animals housed in Micro-Isolator™ Systems, for other housing system applications.



Stack-A-Rack™ Interchangeable See-Through™ Housing Systems meet or exceed ILAR Guide Recommendations for housing rodents.

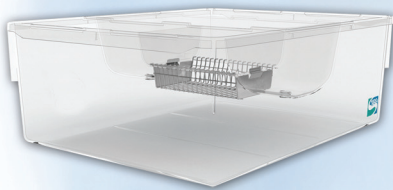

Super Mouse 1800™
(180 in² floor area)



Super Mouse 750™
(75 in² floor area)

- Mount 1 to 8 modular shelves to create your own flexible system
- Rack includes a rigid mobile base, angle uprights and 1 to 8 modular shelves
- Components can be disassembled and cleaned / autoclaved separately
- Two interchangeable cage sizes for housing mice
- Super Mouse 1800™ floor area >180 in² (>1161 cm²)
- Super Mouse 750™ floor area >75 in² (>483 cm²)
- Eliminates wire bar lids for easy observations of animals, food and water
- Watering options: Hydropac™, water bottles or automatic watering

