


# Waste Management Systems

## Mobile Waste Disposal Station

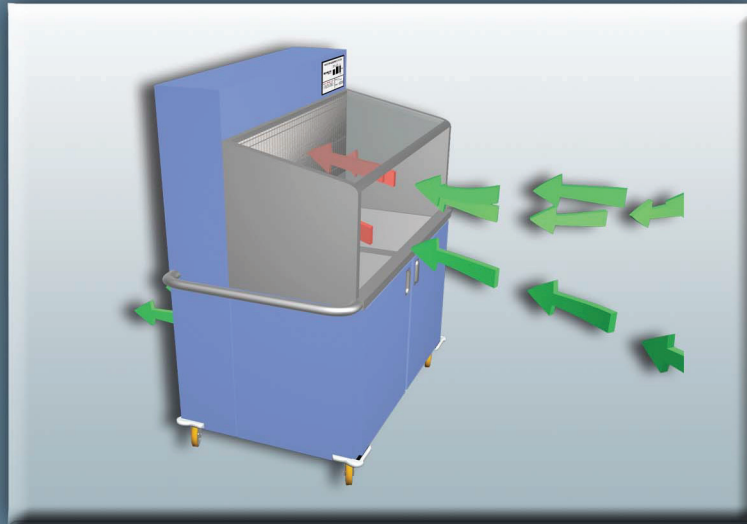

### Manage Waste and Reduce Airborne Particulate.

- Helps Reduce the Spread of Aerosolized Contaminants
- Helps Minimize Personnel Exposure to Particulate
- 3 Stage Air Filtration; Polyester Pre-Filter, Pleated Pre-Filter, HEPA Filter
- HEPA Filtered Air / 99.99% Efficient to .3 Micron Particles
- Access and Change Pre-filters Without Tools
- Work Area Opening (#59020): 37" x 16.50" H (93.98 x 41.91 cm H)
- Dual Stainless Steel Funnel Disposal Ports Channel Waste To Containers
- Holds 2 Plastic Lined Disposal Boxes or Bags (Quick Release Bag Locks)
- Disposal Bags Can be Sealed Within the Unit to Reduce Spread of Particulate
- Pre-Filter and HEPA Filter Warning Lights
- Serves as a Negative Airflow Workstation with Optional Port Cover Plate(s)
- Compact and Mobile Designs for Easy Transport
  
- Overall Dimensions:
  - #59020 46.63" x 31.38" x 66.50" H (118.4 x 79.7 x 168.9 cm H)
  - #59030 52.63" x 35.38" x 70.50" H (133.7 x 89.9 x 179.1 cm H)